




Silicone Ventral Wall Defect Silo Bag

- Designed for staged silo closure of congenital ventral wall defects in infants.
- Silo closed system helps prevent peritoneal fluid loss and bacterial contamination.
- Transparent silicone allows continuous monitoring of bowel condition.
- Open end is supported with a tube of silicone elastomer containing a stainless steel spring.
- Distal end is mesh reinforced.
- Latex-Free

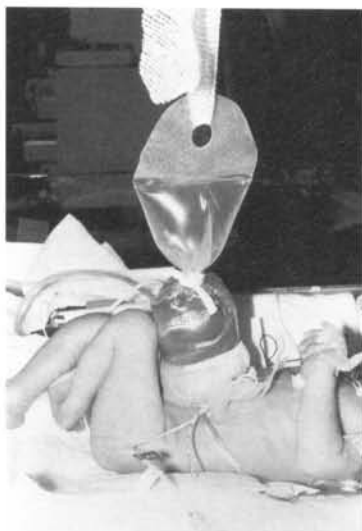

Above: Gastroschisis before Silo insertion.

Right: Gastroschisis with Silo in place.

Photographs courtesy of Loma Linda University Children's Hospital

Silicone Ventral Wall Defect Silo Bag

STERILE

Product Code	Opening Diameter
GR74089-06	3 cm
GR74089-07	4 cm
GR74089-01	5 cm
GR74089-05	6 cm (tapered)
GR74089-02	7.5 cm
GR74089-03	10 cm
GR74089-04	15 cm

Caution: Federal (USA) law restricts these devices to sale or on the order of a physician.
All Bentec Medical, Inc. products are manufactured in the U.S. by Bentec Medical, Inc. in compliance with strict GMP Guidelines.



bentec
MEDICAL

1380 E. Beamer St.
Woodland, CA 95776

Fax: 888-828-8869
Internet: www.bentecmed.com

To order, call 800•767•9175
Visa/Mastercard/Procurement
Cards Accepted